paro[®] sonic

user manual

EN

DE

FR

IT

HE

JA

SWISS
MADE

www.paro.com | info@paro.com

Congratulations on your purchase of the **paro**[®] *sonic* - the new generation of sonic toothbrushes. The hydrodynamic sonic toothbrush with the TOTAL SONIC CARE system - two different toothbrush heads and the interdental function – for thorough and effective tooth care.

The light-weight handle with a vigorous sonic power drive and rechargeable lithium-ion batteries is easy to handle and operate. The timer signals the cleaning process and the travel safety mechanism protects the toothbrush against accidental activation. The **paro**[®] HYDRODYNAMIC SENSITIVE TECHNOLOGY provides for brilliant and gentle cleaning.

Please read this operator's manual carefully before use and keep it for future reference. The operator's manual provides tips and tricks for good and thorough teeth cleaning. Additional information and tutorial videos can be found at www.paro.com under the **paro**[®] *sonic* product.

SAFETY INFORMATION

WARNING

- Please consult your dentist if you have recently undergone an operation or oral treatment.
- Consult your dentist if prolonged or heavy bleeding occurs after using the toothbrush.
- Children or people with intellectual disabilities should be supervised to ensure proper use of the device.
- Ensure that children do not play with the device.

CAUTION

- The power cable cannot be replaced. If the power cable is damaged, dispose of the charger.
- Keep the charger away from water! Never leave the charger above or near a bathtub or washbasin filled with water.
- Connect the charger with the enclosed AC adapter to a wall mounted socket away from the water source and keep away from water.
- Protect the plugged-in charger from hot surfaces.
- The **paro**[®] *sonic* sonic toothbrush complies with the safety standard for electromagnetic devices. However, if you have an implant or pacemaker and as such have questions about the use of the device, please contact your doctor or the manufacturer of the pacemaker or implant.

GOOD TO KNOW


CHOOSE MODE – FOR TOOTHBRUSH HEADS

Press the Start/Stop button (1) again within 2 seconds to switch from SENSITIVE mode to INTENSIVE mode. After 2 seconds, the mode is set and by pressing the button again, the **paro**[®] *sonic* switches off. When switched on again, the **paro**[®] *sonic* starts in the mode that you last used.
-> Memory Effect

Hold the brush head at a 45° angle to the teeth and guide the toothbrush gently from tooth to tooth. Linger for about 2-3 seconds per tooth. Additional wiping movements from the gum to the crown of the tooth (red to white) facilitate the brushing effect.

For each of the four jaw quadrants, guide the toothbrush along the occlusal surface, then along the outer surface, and finally along the inside of the row of teeth (Figure 2).

The quadrant timer assists your brushing and indicates when you should change the quadrant by briefly pausing the toothbrush's function. After four quadrants, the **paro**[®] *sonic* automatically switches off. The toothbrush heads should be changed every 2 months.

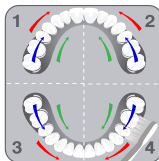① Start/Stop and the mode selection button for the use of the toothbrush heads.

- ① SENSITIVE - 80'000 cleaning movements per minute/2-minute timer
- ① INTENSIVE - 70'000 cleaning movements per minute/3-minute timer

CHOOSE MODE – FOR INTERDENTAL HEADS

With the **paro**[®] *sonic* in INTERDENTAL MODE, you can comfortably clean interdental spaces, braces and implants using the **paro**[®] *isola F* Interdental brushes. Screw the appropriate size of **paro**[®] *isola F* interdental brushes to the interdental head, start the **paro**[®] *sonic* in INTERDENTAL MODE and insert the brush into the interdental space, move it back and forth twice before moving to the next interdental space.

IMPORTANT: Only use the **paro**[®] *isola F* interdental brushes in INTERDENTAL MODE.

② Start/Stop button for the use of the toothbrush heads and travel lock

- ① INTERDENTAL - slight vibrations for cleaning the interdental spaces, implants or braces with **paro**[®] *isola F* interdental brushes

TRAVEL-LOCK

The travel-lock is activated by pressing the INTERDENTAL Start/Stop button for 4 consecutive seconds. This locks the function buttons against a random start.

Pressing the INTERDENTAL Start/Stop button again for 4 consecutive seconds releases the travel-lock. Placing the **paro**[®] *sonic* on the plugged-in charger automatically releases the travel-lock as well.

③ BATTERY STATUS INDICATOR

Before using the **paro**[®] *sonic*, fully charge it for 10 hours. The battery symbol shows you the battery status after brushing your teeth:

Permanent glowing symbol: 100% - 75%

Slowly dimming symbol: 75% - 35%

Rapid blinking symbol: less than 35%

When the symbol is rapidly blinking, the **paro**[®] *sonic* should be charged.

BRUSH HEADS AND INTERDENTAL CLEANING

④ **paro**[®] *sonic* soft-clean - cleans gently, protects sensitive gums and exposed tooth necks

⑤ **paro**[®] *sonic* duo-clean - provides for more intensive cleaning with depth effect

⑥ **paro**[®] *sonic* interdental - with exchangeable **paro**[®] *isola F* interdental brushes, it is the perfect solution to clean niches, interdental spaces and implants as well as braces

⑦ USB charger with brush head holders and power pack

TECHNICAL DATA

Mains voltage: 100-240 V AC, 50 Hz, 3.7V, Class II, IPX7

Battery: Rechargeable lithium-ion battery

Operating time (full to empty battery): ~ 120 minutes

Charging time with first charge: 10 hours

This device complies with EU directives (EMC/RoHS) and low voltage.

DISPOSAL

This device contains a lithium-ion battery. At the end of its service life, the device must be recycled or disposed of in an environmentally responsible manner.

Special containers for disposal and recycling are provided in stores or public collection points.

WARRANTY

paro[®] grants a 2-year warranty from the date of purchase. This warranty applies to any country where this product is distributed by **paro**[®] or a designated distributor.

The following are excluded from the warranty. Damage due to improper use, normal wear and tear, normal use, especially brush heads, and defects that have a negligible effect on the value or operation of the device. The warranty is void if repairs are made by unauthorised persons or the housing is opened. To qualify for the warranty, send the entire device with proof-of-sale device to the official **paro**[®] representative in your country. If this is not known, please contact us directly:

www.paro.com, +41 (0)44 723 11 10 (the usual costs for an international call apply).

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paro[®] *sonic* soft-clean brush head

cleans gently, protects sensitive gums and exposed tooth necks

TOOTHBRUSH MODE - Start/Stop and the mode selection button:

- **SENSITIVE** - 80'000 cleaning movements per minute / 2-minute timer
- **INTENSIVE** - 70'000 cleaning movements per minute / 3-minute timer

INTERDENTAL MODE / TRAVEL-LOCK

INTERDENTAL MODE

For cleaning interdental spaces, braces and implants with **paro**[®] *sonic* interdental head

TRAVEL-LOCK ON / OFF

Press the INTERDENTAL Start/Stop button for 4 consecutive seconds

paro[®] *sonic* duo-clean brush head

provides for more intensive cleaning with depth effect perfect

paro[®] *sonic* interdental head

with exchangeable **paro**[®] *isola F* interdental brushes

paro[®] sonic refill

paro[®] sonic soft-clean #761

paro[®] sonic duo-clean #762

paro[®] sonic interdental #763

paro[®] isola F

#1048

#1041

#1049

#1043

#1042

#1046

#1044

#1047

#1013

Toothpaste recommendation

paro[®] amin #2668

paro[®] sonic friendly & tested

Optimal caries protection
✓ 1250 ppm Amin fluoride

Gum strengthening
✓ SLS free
✓ Panthenol

Mint flavour

paro[®] – developed and produced by

Esro AG | CH-8802 Kilchberg | Switzerland
www.paro.com | info@paro.com